


LSR II Configuration

Hoeschst 33342 Red
670/14


BUV 496
Live/Dead Aqua*

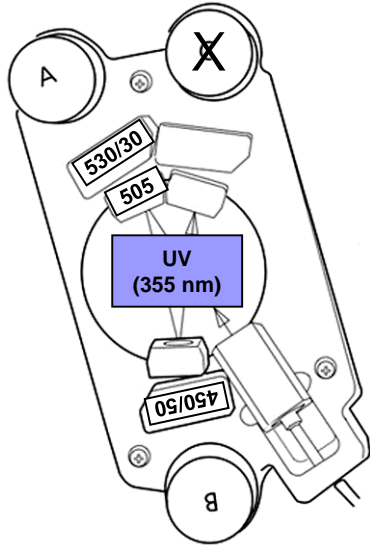

Alexa 350
DAPI*
Hoeschst 33342 Blue
Hoeschst 333258
Live/Dead Blue

BUV395
379/28


APC-H7
APC-Cy7
APC-eFluor 780

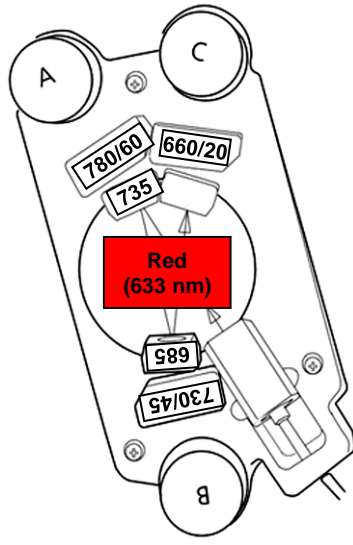

Alexa 700

BV605 610/20 575LP
BV711 710/50 685LP


BV570

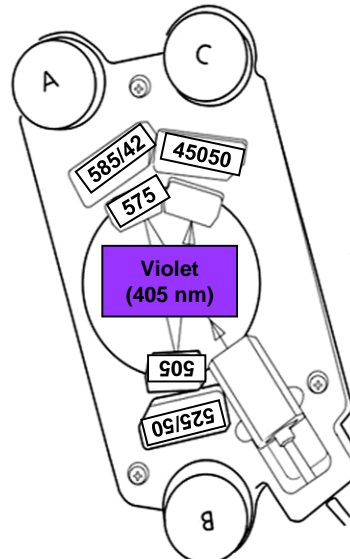

BV 510
BD Horizon V500
AmCyan
CFP
Live/Dead Aqua*

BV570 585/42 575LP
BV605 610/20 575LP

BV 421
Pacific Blue
BD Horizon V450
Live/Dead Violet
DAPI*

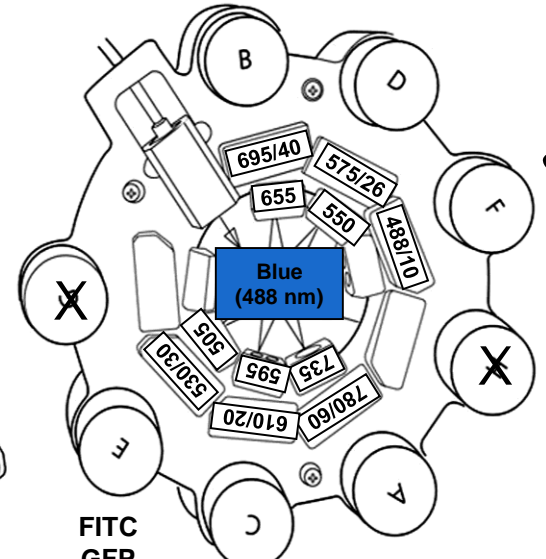

FITC
GFP
Alexa 488
YFP

PI
PE-CF594
PE-Alexa 610
PE-Texas Red

PE-Cy7

7AAD
PerCP
PerCP-Cy5.5
PerCP eFluor 710

PE
Alexa 532
PyroninY
RFP (if not too weak)
tdTomato (if not too weak)

SSC

* Approximate MFI can be detected well by other PMTs